DIODE LASER LX16 PLUS

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
- TMJ TREATMENT

FUNCTION TO USE	 PAIN THERAPY
TIP TO USE	 DT-30
POWER	4 W

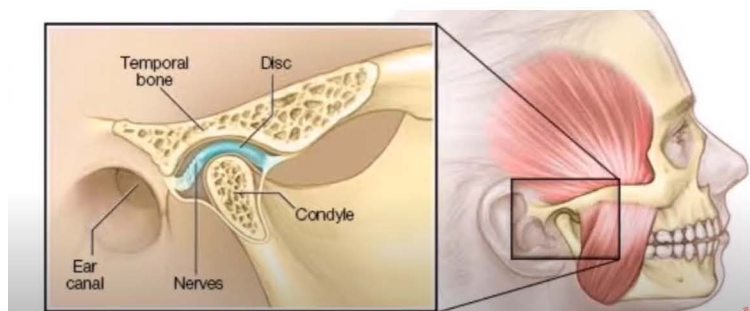
- ALL TIPS USAGE ON FUNCTIONS

HEAD TIPS	FUNCTIONS TO USE WITH TIP	BEST FOR
 DT-15	▪ Aphthous Ulcers: Effective for speeding up healing and reducing discomfort associated with canker sores. ▪ Herpes: Helps reduce inflammation and pain in herpes lesions in and around the mouth. ▪ Pain Therapy: Ideal for relieving pain, especially in areas like sensitive teeth, gums, or after dental procedures. ▪ Dentin Hypersensitivity: Effective for reducing sensitivity, especially from hot, cold, or sweet stimuli. ▪ Burning Mouth Syndrome: Relieves burning sensation in the mouth.	▪ Oral pain management, ulcers, herpes, burning mouth syndrome, and dentin hypersensitivity.
 DT-20	▪ Wound Healing: Encourages faster healing of oral wounds, ulcers, and small post-surgical areas. ▪ Pain Therapy: Can also be used to relieve pain associated with dental procedures, TMJ pain, or post-treatment discomfort.	▪ wound healing, and pain relief (especially in deeper tissue or larger areas).
 DT-30	▪ Haemostasis: This head is primarily used for controlling bleeding after minor dental procedures or injuries. ▪ Wound Healing: It can also help promote faster healing of oral wounds or post-surgical tissue repair. ▪ Pain Therapy: May be used for deeper tissue pain, particularly around the gums or jaw area.	▪ Stopping bleeding after procedures, wound healing, and deep tissue pain relief. Not used for teeth whitening.
 BT-8 TIP	▪ Teeth Whitening: One of the most effective tips for teeth whitening when paired with whitening gels. It helps speed up the whitening process by directing energy to the teeth's surface. ▪ Dentin Hypersensitivity: Can be used to treat teeth that are sensitive to hot, cold, or sugary foods by desensitizing the nerves. ▪ Pain Therapy: Also useful for mild pain relief, particularly in cases of tooth sensitivity or gum discomfort.	▪ Teeth whitening, dentin hypersensitivity, and mild pain relief.
 WHITENING-4 TIP	▪ Teeth Whitening: Designed specifically for whitening treatments. ▪ Desensitization: Helps reduce sensitivity during and after whitening treatments.	▪ Teeth whitening and desensitizing teeth during or after whitening procedures.

POWER USAGE DIFFERENT

- 1. PULSED MODE (E.G. 5 HZ, 65% DUTY):**
 - THE BEAM IS ON AND OFF REPEATEDLY.
 - EACH PULSE ENDS, THEN A SHORT REST PERIOD HAPPENS MEANING THE LASER DOESN'T STAY FIRING CONTINUOUSLY AT ONE SPOT.
 - THIS GIVES TISSUE A CHANCE TO COOL DOWN AND PREVENTS EXCESSIVE CARBONIZATION.
- 2. CW MODE (CONTINUOUS WAVE):**
 - THE BEAM STAYS ON CONTINUOUSLY FOR THE FULL SET TIME.
 - NO OFF BREAKS UNTIL THE TIMER RUNS OUT (OR YOU STOP IT).
 - SAFE ONLY AT VERY LOW POWER (MILLIWATTS), USED FOR THERAPY/BIOSTIMULATION.

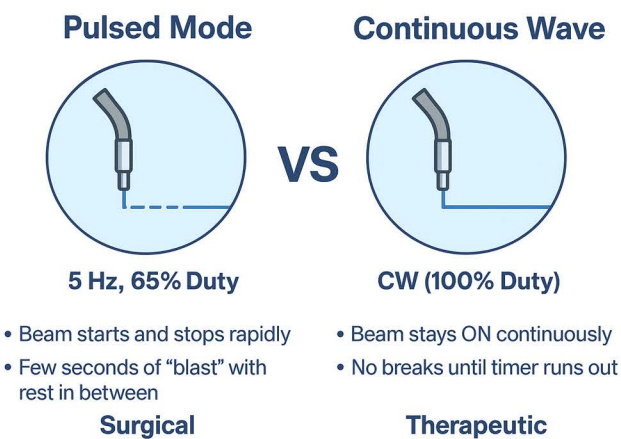SPECIFICATION

WAVELENGTH	POWER (P _{MAX})	CLASS	MODE	MAIN USE
450 ± 20 NM (BLUE)	3 W	CLASS IV	CW / PULSED (1 HZ–20 KHZ)	DISINFECTION, BIOFILM REMOVAL, SUPERFICIAL COAGULATION
650 ± 20 NM (RED)	200 mW	CLASS II	CW	PHOTOBIMODULATION (HEALING, PAIN RELIEF, NERVE THERAPY)
976 ± 20 NM (INFRARED)	5 W	CLASS IV	CW / PULSED (1 HZ–20 KHZ)	SOFT TISSUE SURGERY (PULPOTOMY, FRENECTOMY, COAGULATION, STERILIZATION)
650 NM (AIMING BEAM)	< 5 mW	CLASS II	CW	VISUAL GUIDE ONLY
BATTERY	11.1 V, 2600 MAH ×2 (57.7 WH)			PORTABLE OPERATION